



USER MANUAL



Pedal Exerciser (EG0087)

Copyright Zoom Mobility

User Manual Booklet 2023

Pictures used are for illustration purposes only

Please read the instructions carefully before using this product.

The Zoom Mobility Pedal Exerciser is a compact exercise bike that stores conveniently out of the way when not in use. Use it regularly to stimulate blood circulation and increase muscle strength: it is a great way to enjoy moderate cardiovascular workout. Place it on a tabletop to exercise your arms or on the floor to enjoy exercising your legs from your chair.

WARNING! Do not stand on your pedal exerciser.



Preface

- I. Open the package and take out the product and all components.
- II. Check each part carefully before use. If anything becomes broken, do not use.
- III. Keep these instructions in a safe place for future reference.

How To Assemble


Examine the parts included in the packaging to make sure they are all there.

1. Pull firmly on the ring at the joint of the frame leg.
2. Pull the leg out until it click into place, repeat with the other leg.

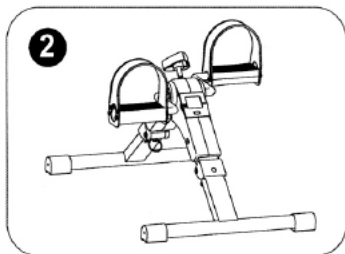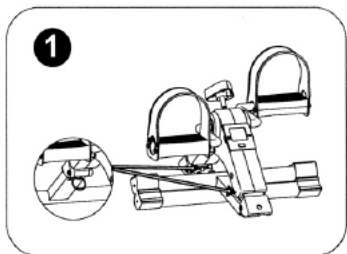


Monitor Functions


1. Loosen and remove the screws on each side of the screen unit, then carefully remove the screen unit from its socket.

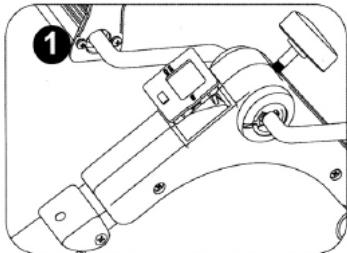


2. On the back of the screen unit, lift up the battery tabs and side out the battery.

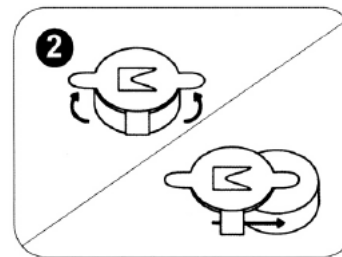

3. Insert a new 1.5VAG13 battery positive side down and close the tabs.

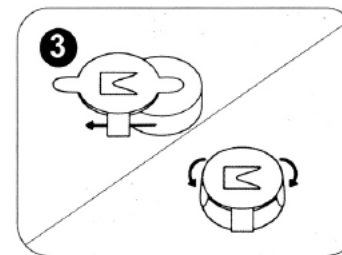


4. Replace the screen unit in its socket.

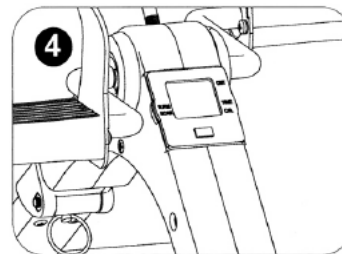

Display Warnings


To turn on the display, press once on the red button or start pedaling. The display automatically shuts off after 4 minutes of inactivity. All the display functions are activated as soon as you begin pedaling and pause when you stop. To toggle between them, press the red button until the arrow points to the desired function.

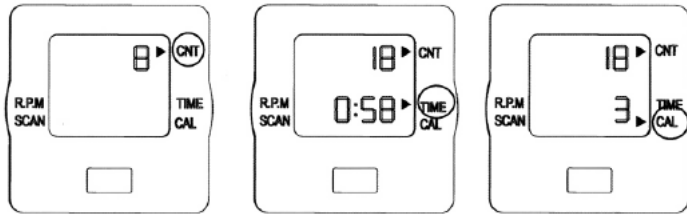

CNT:
Number of rotations completed.

TIME:
Workout Duration.

CAL:
Estimated calories burned.

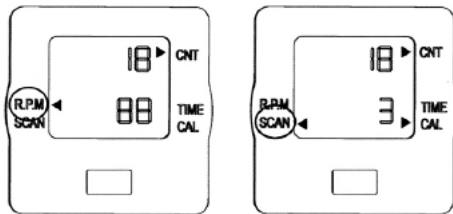

R.P.M.:
Rotations per minute.

SCAN:
Toggle (TIME > CNT > RPM > CAL).

The guarantee does not extend to those items which may need replacement due to normal wear and tear, or damage to the product caused by misuse or accident for which Zoom Mobility or its agent cannot be held responsible. This warranty does not include service calls.

Warranty

One-Year Warranty

- Pedal Exerciser (Digital Display)

Your Zoom Mobility product is fully guaranteed for one (1) year from the date of purchase against faults arising due to defects in the manufacture or materials. This warranty does not detract from it but is in addition to your legal rights.

Any attempt to open or dismantle this item renders the guarantee void on that item.

Not Covered Under Warranty